



FOR EXPORT ONLY

The Advancement of **NEW**
RF Technology

Advanced Anti-aging RF Technology

SECRET is an intralesional Fractional Radio Frequency (RF) device, that utilises minimally invasive micro-needles for non surgical face lifting, wrinkle reduction and scars for dermal coagulation.

Compared with all fractional laser systems, Secret can treat all dark skin types with no risk of skin burns and PIH making the treatment faster with relatively no down time.

Secret^{RF}

Micro-Needle Fractional RF System

Indications

Face Treatment

- Non-surgical Face Lifting
- Skin Tightening
- Pore Reduction
- Wrinkle Reduction
- Skin Rejuvenation (Whitening)
- Acne Scars

Body Treatment

- Scars
- Stretch Marks



Secret^{RF}

Micro-Needle Fractional RF System

⦿ No down time, Painless

Precisely controlled RF energy is delivered only to target areas by selectively adjusting the depth of the needles from 0.5mm to 3.5mm, with an advanced penetration mechanism for dermal coagulation.

⦿ Optimal and uniform thermolysis

Ideal distance among needles allows applying uniform bipolar RF energy in the target areas to coagulate evenly without any risk of energy interference caused by overlapping of RF energy in the close distance.

⦿ Precise and safe treatment

Light variation function depending on each operation step allows practitioners to track overall treatment procedure during treatment for patient safety as well as better result.

⦿ Fast and convenient operation

10.4" full color touch interface and various modes for inspection and cleaning are user friendly and easy to operate.

Specification


Micro-Needle Fractional RF System	
Max. Power	20W ± 20%
RF Frequency	2MHz ± 20%
Intensity	0 ~ 100 Level (10step)
Repetition	0.2 / 0.5 / 1 / 2sec / Single
Needle Thickness	0.3mm
Needle Figure	25 electrode needles (5×5), 7.8mm×7.8mm / pitch 1.95mm (Standard) 16 electrode needles (4×4), 7.8mm×7.8mm / pitch 2.6mm (Optional)
Treatment Area	1cm×1cm
Adjustable Depth	0.5 ~ 3.5mm (0.1step)
Treatment Duration (RF emission)	50ms ~ 950ms
Electrical Power	100 ~ 240VAC, 50 / 60Hz, 220VA under
MONO	
Max. Power	50W ± 20%
RF Frequency	2MHz ± 20%
RF Electrode	Mono polar
Repetition	CW
RF Electrode Size	φ15



Sarkis & Bakhos Strt, Sector 2
GETO GROUP Bldg, 4th Floor
Jdeidet El-Metn - Beirut, Lebanon
Tel + 961 1 876 579
Fax + 961 1 876 580
Mobile + 961 3 945 738


www.medistatintl.com

How Secret works?

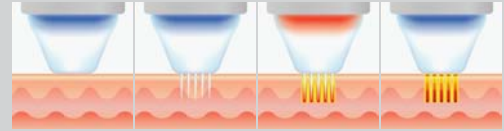


Micro-needle Cartridge is placed on the skin surface. Micro-needles penetrate inside of the skin. Bipolar RF energy is directly applied to the target areas. Collagen regeneration and new elastic production process begins.

Benefits of Micro-needle Fractional RF

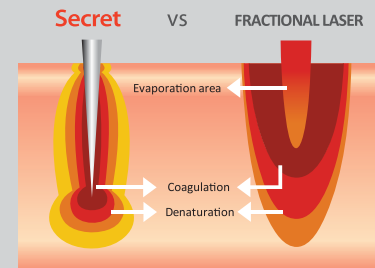

Ideal distance for uniform thermolysis

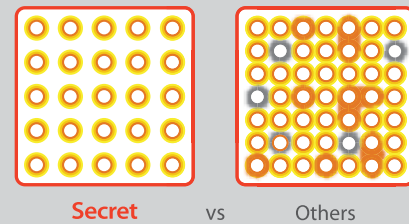

Light variation for precise treatment



Mono Polar RF handle

