



Can operate
2 main handpieces
at once

Optional Large Hand-piece


Applicators	<ul style="list-style-type: none"> ● 2 medium size handpieces 108mm(W) x 188mm(D) x 138mm(H) * Small handpiece (Optional) * Large size hand piece (Optional) 118mm(W) x 220mm(D) x 142mm(H)
Pulse type	<ul style="list-style-type: none"> ● Pulsed & Continuous
Cooling Temperature	<ul style="list-style-type: none"> ● -10°C ~ 5°C
Treatment time	<ul style="list-style-type: none"> ● 01:00 ~ 90:00 min
Vacuum Level	<ul style="list-style-type: none"> ● Level 1 to Level 5
Input Power	<ul style="list-style-type: none"> ● 100~110V, 220 ~ 240V,50/60Hz (max 600w)
Panel	<ul style="list-style-type: none"> ● 10.2 inch LCD & Touch screen
Dimension (WxDxH mm)	<ul style="list-style-type: none"> ● 1,066x379x438
Weight	<ul style="list-style-type: none"> ● 60 Kg

Medistat international s.a.r.l
Sarkis & Bakhos street, sector 2, Geto Group Bldg,
4th Floor, Jdeidet el - Metn, Beirut, Lebanon,
Tel: +961 | 876579 Fax: +961 | 876580

info@medistatintl.com
www.medistatintl.com



COOL
SHAPING™

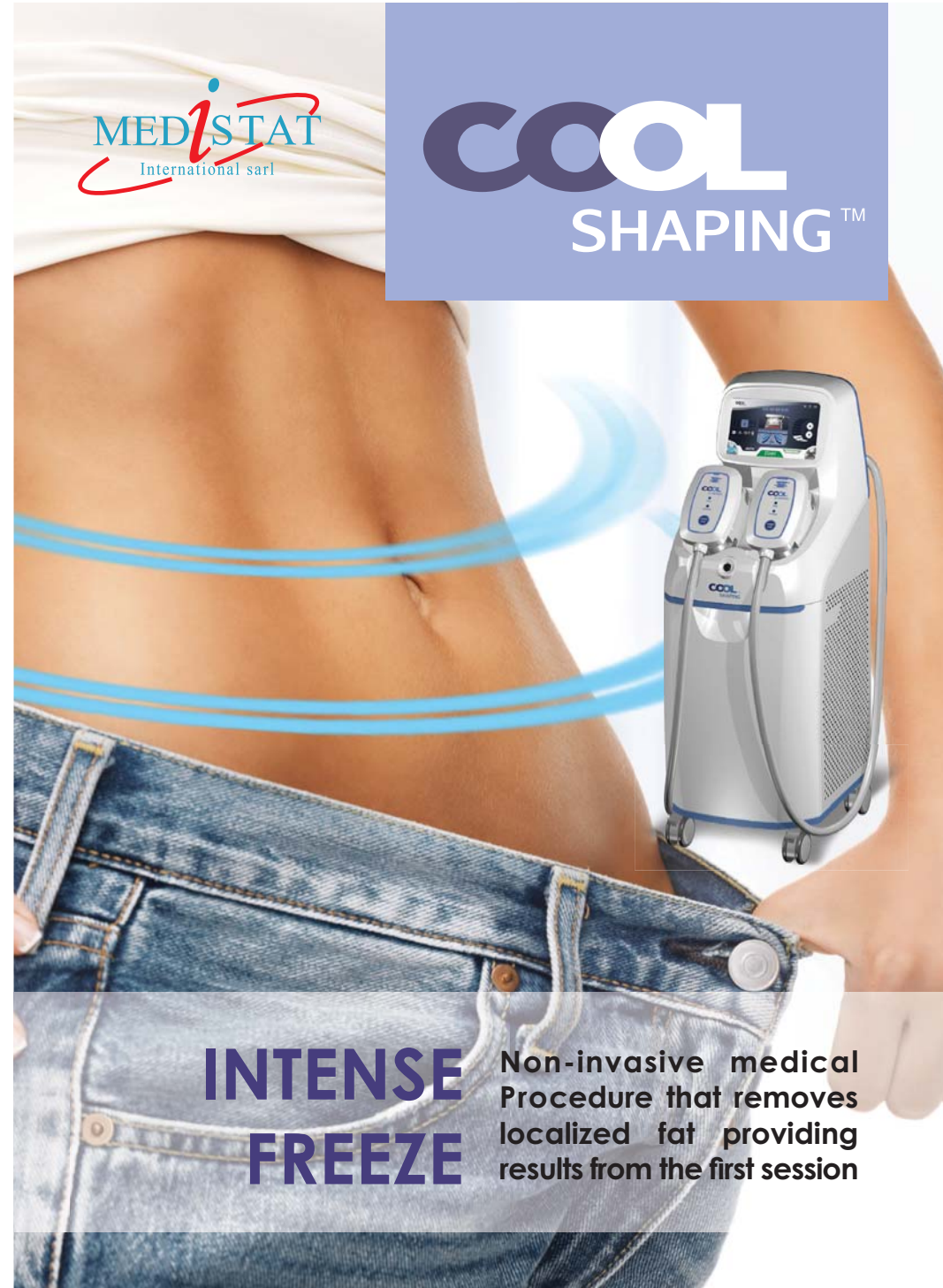

**INTENSE
FREEZE**

Non-invasive medical
Procedure that removes
localized fat providing
results from the first session




Cool Shaping™ offers the most advanced body slimming solution by the integrated technology™ of cryo, vacuum and led.

GET BETTER RESULTS FOR BODY SLIMMING TREATMENT by cryolipolysis technology

Integrated Technology™ by Cryo, Vacuum and LED

How CoolShaping™ works?

CoolShaping employs selective cryolipolysis (Fat freezing) using a set of plates to reach temperature up to -10°Celsius. An anti-freeze pad is placed over the skin fold so that the skin exposed to the freezing plates will not be damaged. The fat cells disappear after apoptosis (programmed cell death) and are removed completely by macrophages in a slow, gradual and progressive way.

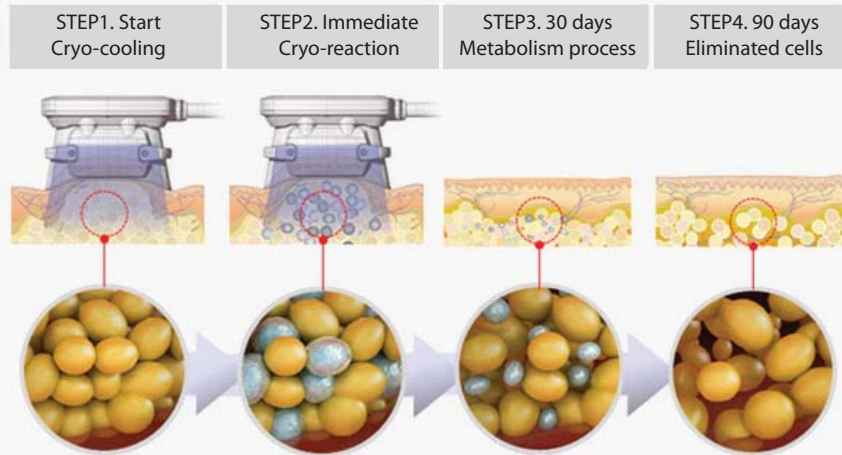

Apply Coolshaping applicator at the target area. Start cryo-cooling to reach the needed temperature.

The fat cells in the target area begin responding to cryo-crystallization

Following the treatment fat cells reach a biological death and removed by metabolism process


The fat cells are eliminated

CRYO

VACUUM

LED

Cool Shaping™ offers the most advanced body slimming solution by an integrated technologies™ cryo, Vacuum and Led. The system delivers cryo-energy during the treatment with vacuum and LED technologies which enhances the treatment effects.

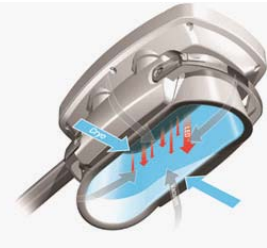

STRONG CRYOLIPOLYSIS SYSTEM

- Two powerful cooling Peltiers attached to the handpiece.
- Automatically control the cooling temperature.
- Maximum -10 °C



630nm Red LED



STABLE VACUUM TECHNOLOGY

- Controllable vacuum level: step 1 (low) to step 5 (high)
- Various vacuum pulses: Pulsed and continuous mode

BEFORE

AFTER



BEFORE



AFTER

